



**SANARE
MEDICAL**

Innovative technologies, products and
services to improve human life.

FibroScan® Mini+ 430

The mobile non-invasive solution for liver disease management

Designed for multi-site sharing and near-patient testing. Powered by unique, patented and validated parameters LSM by VCTE™ and CAP™ for liver fibrosis and liver steatosis assessment.

Ideal point-of-care solution for large sized and community hospital-based settings.



NEW Guided VCTE™

The next generation VCTE™ technology

Intuitive features for faster examinations and simplified scanning¹

- Two new indicators that allow for quick identification of the optimal measurement location
- AutoScan feature which triggers 10 valid measurements through a single click
- SmartExam features fully included: SmartDepth and Continuous CAP™



Stiffness indicator



CAP™ indicator

Smartly adaptable to all settings

- Light and easy to handle (11 lbs./5 kg) with 12.1 inch touch screen
- Battery-powered device
- Two probe connectors to easily switch between probes during exam

Adapted to all patient morphologies

Probes S+, M+ and XL+



Transport suitcase

Airplane-compatible



Enhance your FibroScan® experience with the **Liver Health Management platform**: an innovative cloud-based solution to streamline your workflow



Winner of the **Red Dot Design Award** for innovative and ergonomic design

What makes FibroScan® unique?

- A painless exam that can be performed in **less than 4 minutes** with immediate results at the point-of-care.¹
- Can be performed by **any trained operator** (physician, nurse, or medical assistant).
- Standardized examination with exceptional precision and reproducibility that can be utilized in **99% of patients**.²
- An **ecosystem of solutions** developed by Echosens to support clinical decisions for physicians: Liver Health Management platform, Scores by Echosens^{3,4}, Interpretation Guide, myFibroScan, FibroScan® Gateway and educational support.



Renowned Publication Presence & Endorsement in Clinical Practice Guidelines

4,200+ peer-reviewed publications

180+ international guidelines

More than 4,200 peer-reviewed publications and 180 international guidelines advocate the use of FibroScan® as the reference non-invasive solution for liver fibrosis, cirrhosis, and steatosis assessment across all etiologies of chronic liver disease (viral hepatitis, MASLD/MASH, alcoholic liver disease).^{5,6,7}

AASLD

ADA

AHA

APASL

EASL

AACE


AGA

AISF

Baveno VII

NICE Guidance

Examinations with FibroScan® can inform treatment decisions across the spectrum of disease

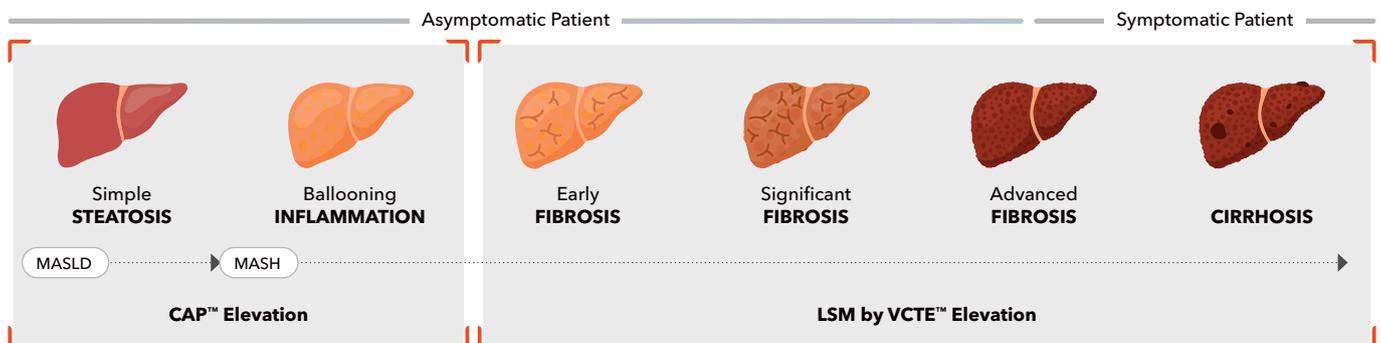

Interested in FibroScan® Mini+ 430 for your practice?
Contact us on info@sanaremedical.com

LSM: Liver Stiffness Measurement / VCTE™: Vibration Controlled Transient Elastography / CAP™: Controlled Attenuation Parameter / MASLD: Metabolic dysfunction-associated steatotic liver disease (formerly known as NAFLD) / MASH: Metabolic dysfunction-associated steatohepatitis (formerly known as NASH) / LHM: Liver Health Management

References

1. Based on internal data. E431R020.1. 2. Myers, Robert P et al. "Feasibility and diagnostic performance of the FibroScan XL probe for liver stiffness measurement in overweight and obese patients." *Hepatology* (Baltimore, Md.) vol. 55,1 (2012): 199-208. doi:10.1002/hep.24624. 3. Newsome, Philip N et al. "FibroScan-AST (FAST) score for the non-invasive identification of patients with non-alcoholic steatohepatitis with significant activity and fibrosis: a prospective derivation and global validation study." *The Lancet. Gastroenterology & hepatology* vol. 5,4 (2020): 362-373. doi:10.1016/S2468-1253(19)30383. 4. Sanyal, Arun J et al. "Enhanced diagnosis of advanced fibrosis and cirrhosis in individuals with NAFLD using FibroScan-based Agile scores." *Journal of hepatology* vol. 78,2 (2023): 247-259. doi:10.1016/j.jhep.2022.10.034. 5. Baveno VII - R. de Franchis et al. "Renewing consensus in portal hypertension." *Journal of hepatology* vol. 76,4 (2022): 959-974. doi:10.1016/j.jhep.2021.12.022. 6. European Association for the Study of the Liver. "Electronic address: easloffice@easloffice.eu et al. "EASL Clinical Practice Guidelines on non-invasive tests for evaluation of liver disease severity and prognosis - 2021 update." *Journal of hepatology* vol. 75,3 (2021): 659-689. 7. Kanwal, Fasiha et al. "Clinical Care Pathway for the Risk Stratification and Management of Patients With Nonalcoholic Fatty Liver Disease." *Gastroenterology* vol. 161,5 (2021): 1657-1669.

The FibroScan® 430 Mini+ is intended to measure liver stiffness (E) using Vibration Controlled Transient Elastography (VCTE™) at 50 Hz shear wave frequency and liver ultrasound attenuation coefficient (CAP™)* at 3.5 MHz. FibroScan liver stiffness measurements (LSM) by VCTE™ may aid the physician in determining the likelihood of cirrhosis and may be used, taken in context with other clinical and laboratory data, as an aid in the assessment of liver fibrosis. FibroScan CAP™ measurements may be used, taken in context with other clinical and laboratory data, as an aid in the assessment of hepatic steatosis. FibroScan® is indicated as a non-invasive aid to clinical management, diagnosis, and monitoring of adult and pediatric patients with confirmed or suspected liver disease, as part of an overall assessment of the liver. Results in the pediatric population should be interpreted while considering the clinical condition and the overall patient profile. The FibroScan® device is intended for use by healthcare professionals in hospitals, clinics or any facility where healthcare is provided. * CAP™ refers to ultrasound attenuation coefficient (originally defined as Controlled Attenuation Parameter). CAP™ on S+ probe is only available with SmartExam capability. © 2024 Echosens - Echosens™ and FibroScan® are trademarks owned by Echosens SA. All rights reserved. One-pager FibroScan® Mini+ 430 - v1 USA - 12/2023.





**SANARE
MEDICAL**

SANARE MEDICAL PVT LTD

Website : www.sanaremedical.com

Email : info@sanaremedical.com

Tel : +91 9810275300 , +91 98101 37486

Corp. Office: E-354 Greater kailash Part 1
New Delhi 110048

Regd. Office : DP 61 Phase 1 Udyog Vihar
Gurugram Haryana 122001